


Transfer Solution to Qubit Tubes

- A1: 1000ul tiprack-our tiprack needs to sit on clear base to fit in slot
- A2: Qubit tube container (WITH black spacer)- leave first well blank then skip a well for every tube placed. Skip a row in between wells to leave a space for caps.
- A3: Trash container
- B1: 15ml round bottom tube containing qubit solution (custom container)
- B2: Overflow tube container for runs greater than 24 samples (Use black spacer)

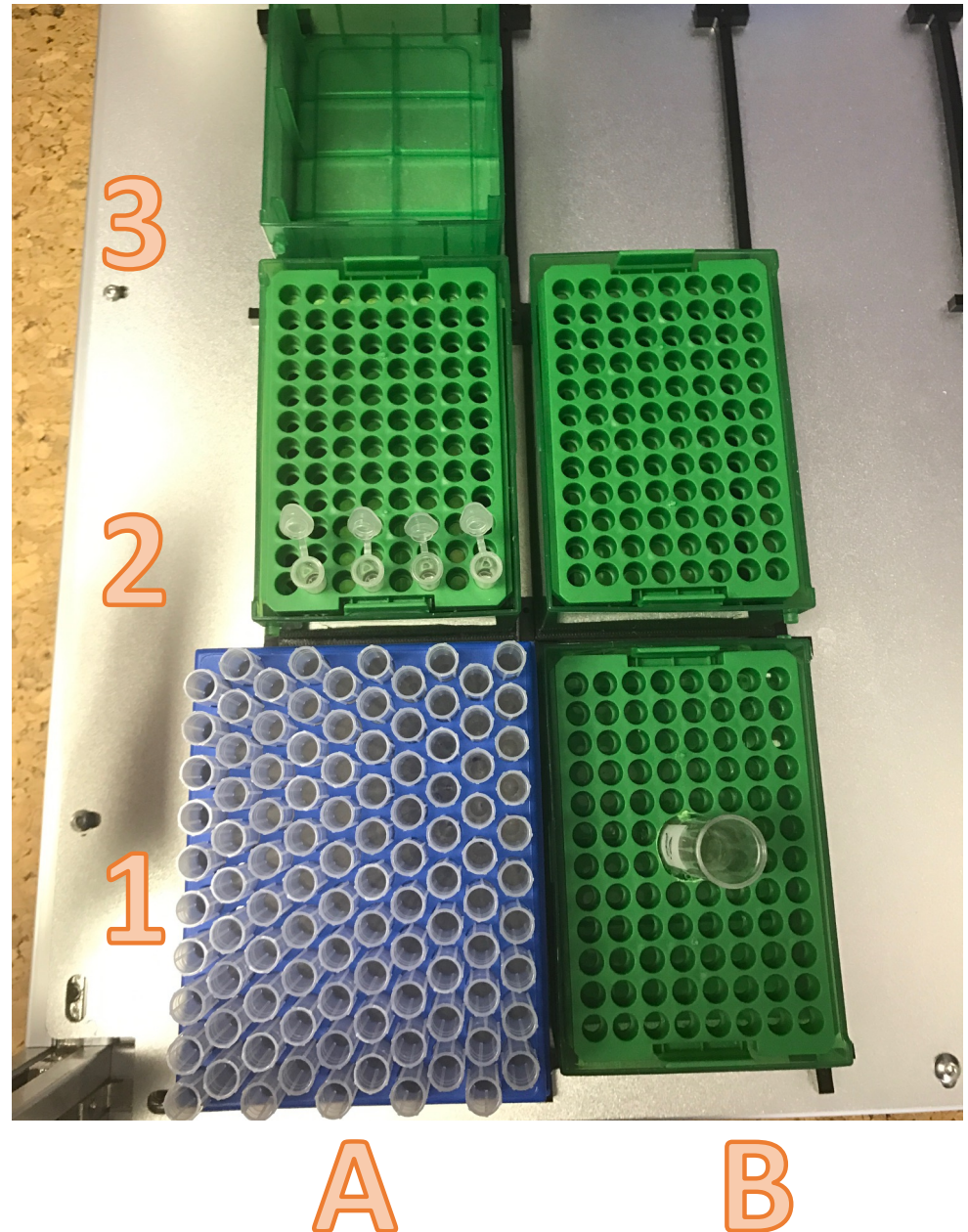

Transfer Samples to Qubit Tubes

D1: Sample Container- Samples in PCR strips- skip row between samples to leave room for caps

D2: Qubit tube container (WITH black spacer)- leave first well blank then skip a well for every tube placed. Skip a row in between wells to leave a space for caps. Has qubit solution from previous step.

D3: Trash Container

E1: 10ul tiprack- take out the tips in first, third, fifth, and seventh column

E2: Overflow tube container for runs greater than 24 samples

E3: If standards needed, place two 0.5ml tubes containing the standards in slot B1 and D1 of green container

