


GLP documentation

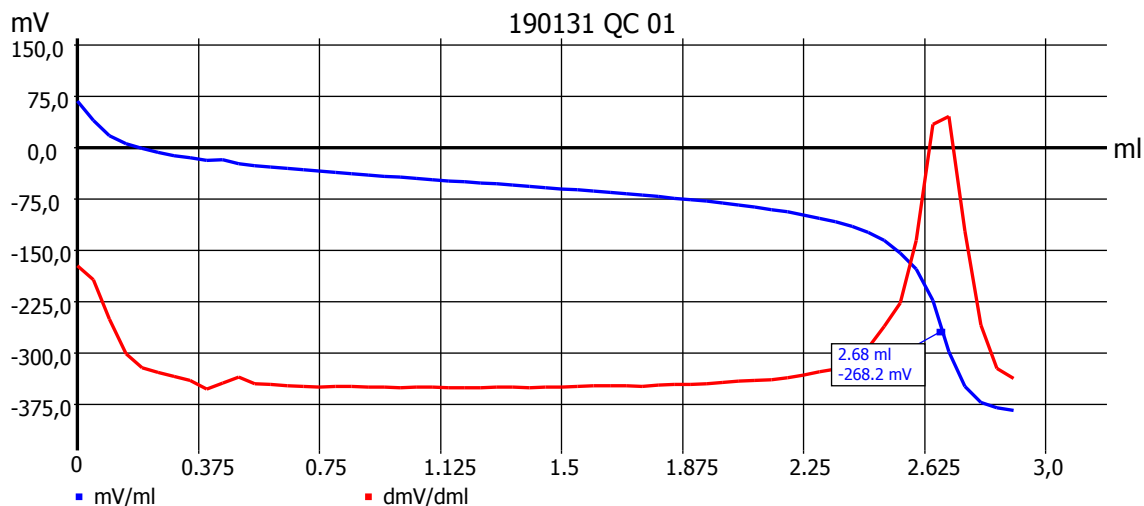

Method data

Method name:	TAN ASTM 664	Titration duration:	11 m 59 s
End date:	31.01.19	End time:	10:58:53

Titration data

Sample ID:	190131 QC 01	Weight:	25.36000 g
Start mV:	68.5 mV	End mV:	-383.6 mV

EQ:	2.678 ml / -268.2 mV	TAN mg KOH/g:	0.555
-----	----------------------	---------------	-------

Calculation formula

TAN mg KOH/g:	$(EQ1-B)*T*M*F1/(W*F2)$	Mol (M):	56.10000
Blank value (B):	0.0000 ml	Titre (T):	0.0937 (a)
Factor 1 (F1):	1.0000	Weight (W):	25.36000 g (m)
Factor 2 (F2):	1.0000	Statistics:	Off

Device information

Device: TitroLine 7750
 Serial number: 10018800
 Software version: 1706b

Method data overall view

Method name:	TAN ASTM 664	Created at:	11/29/17 17:11:00
Method type:	Automatic titration	Last modification:	11/29/17 17:11:00
Measured value:	mV	Damping settings:	strong
Titration mode:	Linear	Documentation:	GLP
Linear steps:	0.050 ml		
Measuring speed / drift:	User-defined:	minimum holding time:	07 s
		maximum holding time:	20 s
		Measuring time:	04 s
		Drift:	10 mV/min
Initial waiting time:	10 s	Stirrer control:	free
Titration direction:	Decrease		
Pretitration:	Off		
End value:	Off		
EQ:	On (1)		
Slope value:	Flat	Value:	120

Dosing parameter

Dosing speed:	100.00 % (20.00 ml/min)	Filling speed:	30 s
Maximum dosing volume:	6.00 ml		

Unit values

Unit size:	10ml
Unit ID:	10036631
Reagent:	ming
Batch ID:	no entry
Concentration [mol/l]:	0.09370
Determined at:	01/31/19 10:23:41
Expire date:	--
Opened/compounded:	01/01/00
Test according ISO 8655:	01/01/00
Last modification:	01/31/19 10:23:46

Device information

Device: TitroLine 7750
Serial number: 10018800
Software version: 1706b