

Run Info

Experiment Name	200123-CLIMARES
Sample ID	MTPX6-6-M40-S00-L20
Run ID	667fc10b-0965-438a-be35-558267cf6e55
Flow Cell Id	PAE48907
Start Time	January 23, 16:01
Run Length	3d 0h 5m

Run Summary

Reads Generated	5.6 M
Bases Generated	73.64 Gb
Estimated Bases	77.77 Gb


Run Parameters

Flow Cell Type	FLO-PRO002
Kit	SQK-LSK109
Basecalling	on
Specified Run Length	72 hours
Initial Bias Voltage	-165 mV
FAST5 Output	Enabled
FAST5 Output Options	zlib_compress,fastq,raw
FAST5 Reads per File	4000
FASTQ Output	Enabled
FASTQ Reads per File	4000
Active Channel Selection	Enabled
Mux Scan Period	1 hour 30 minutes
Reserved Pores	0 %
Basecall Model	High-accuracy basecalling
Barcode Kits	EXP-NBD104
Trim Barcodes	off


Versions

MinKNOW Core	3.6.1
Bream	4.3.16
Guppy	3.2.8

Cumulative Output Reads

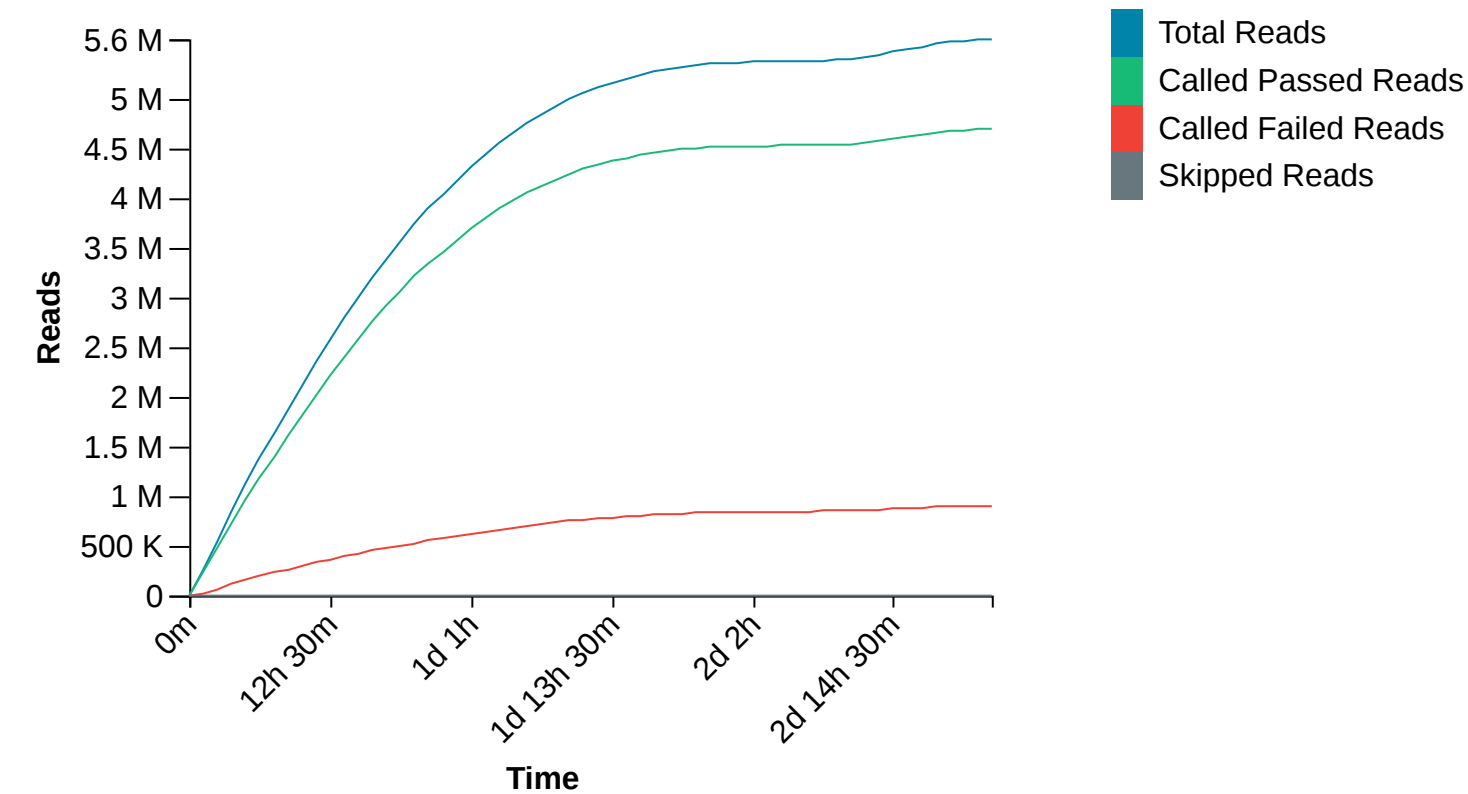


Cumulative Output Bases

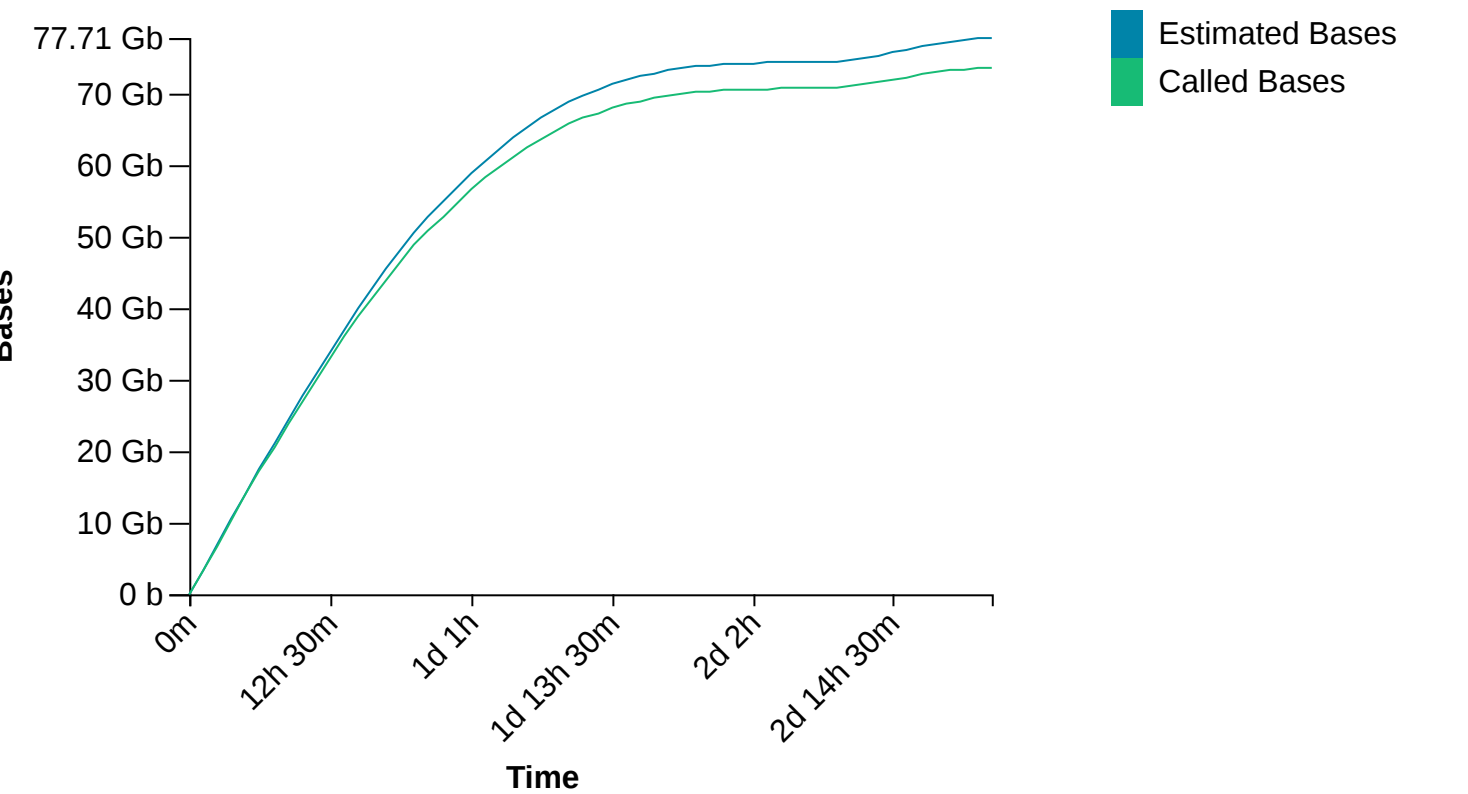


Read Length Histogram Estimated Bases

Estimated N50: 24.35 Kb

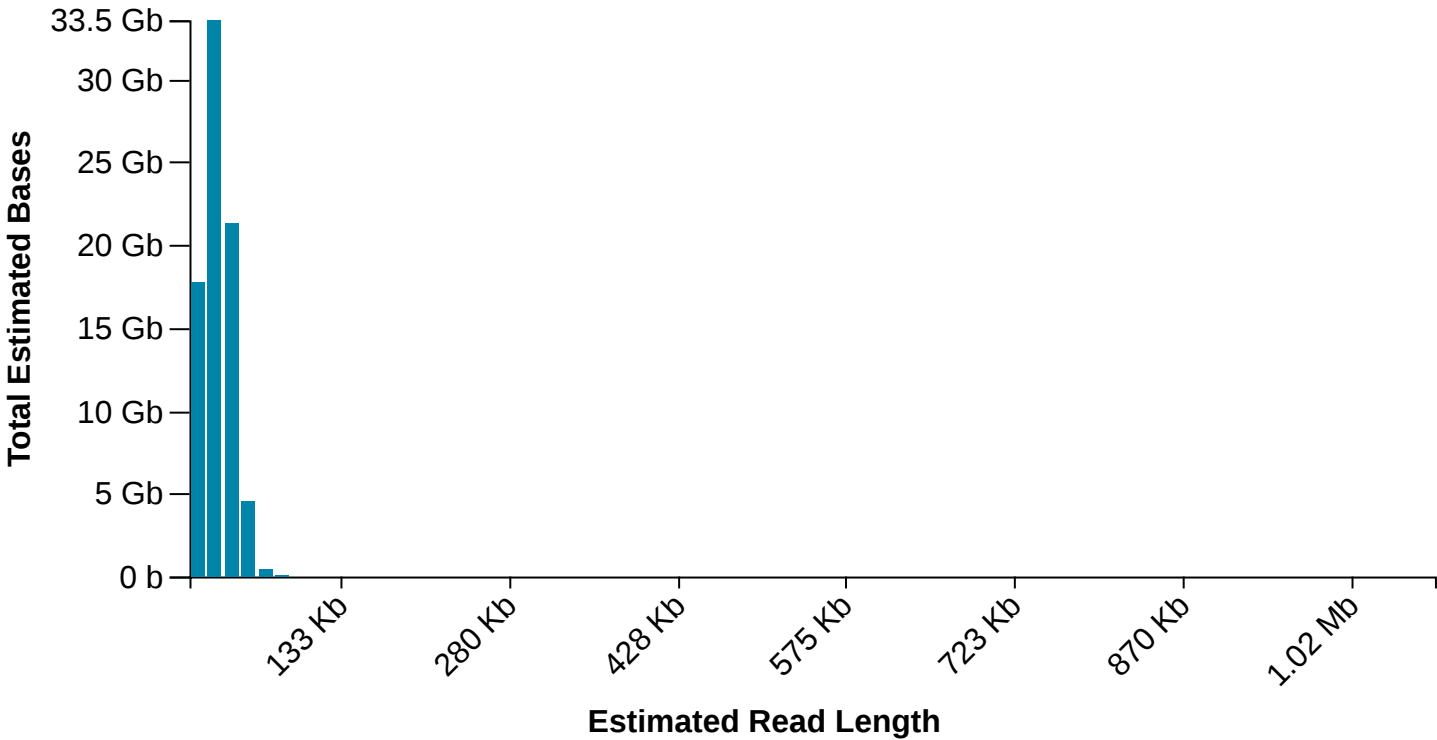


Read Length Histogram Basecalled Bases


Estimated N50: 23.51 Kb

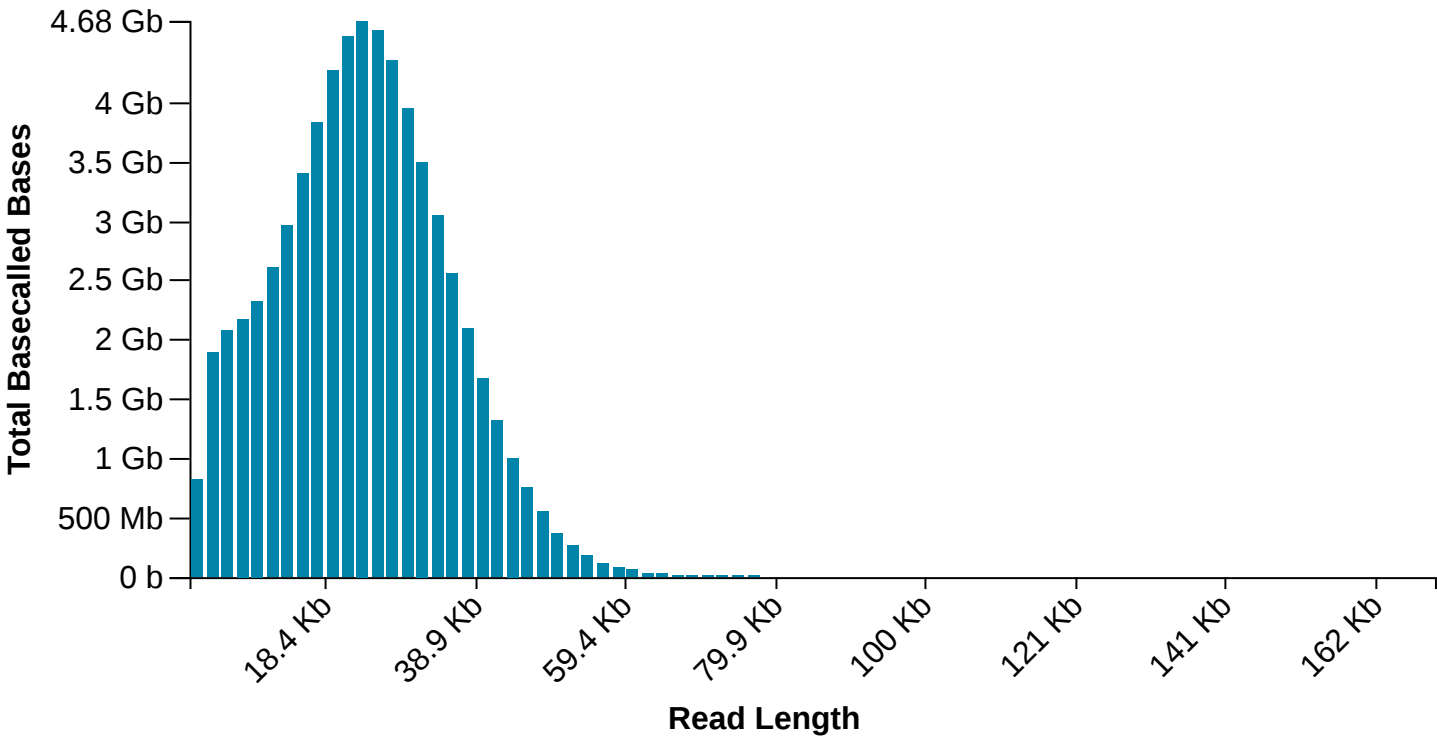


Barcode Read Counts

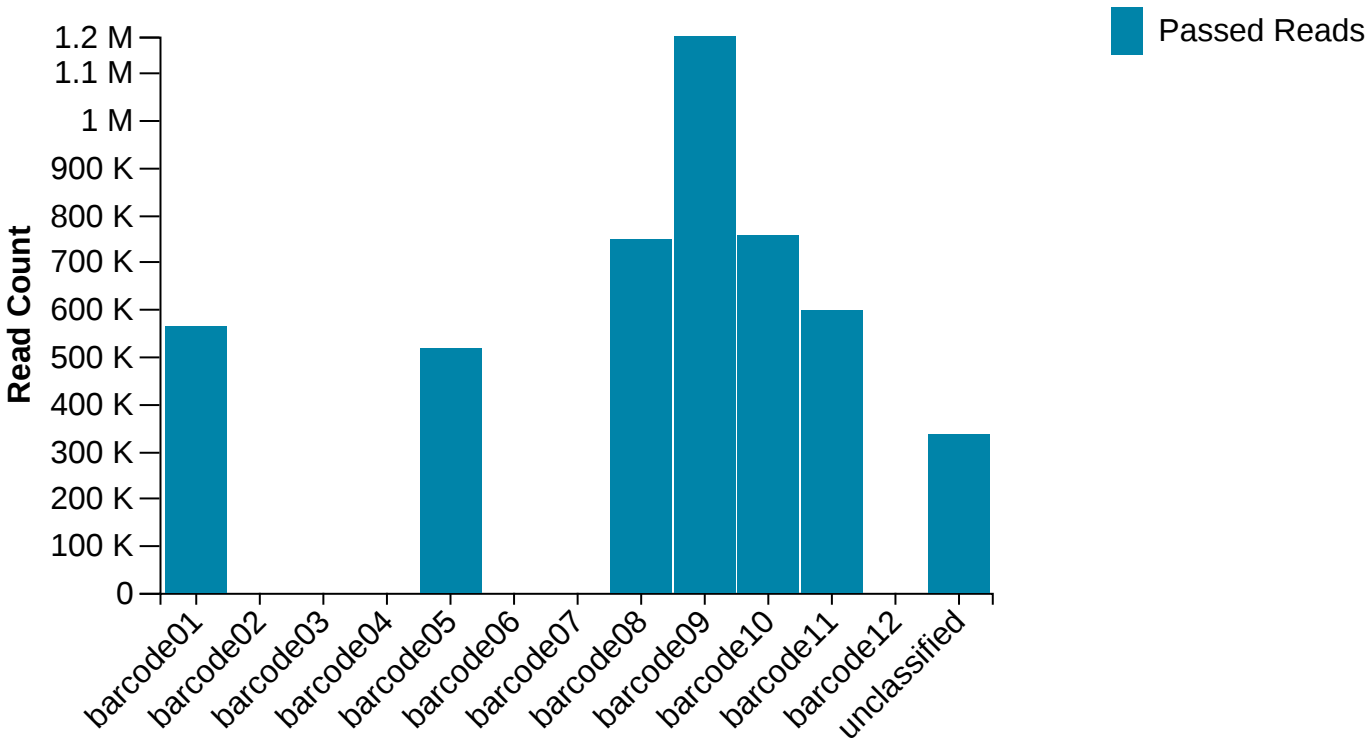


Duty Time Grouped

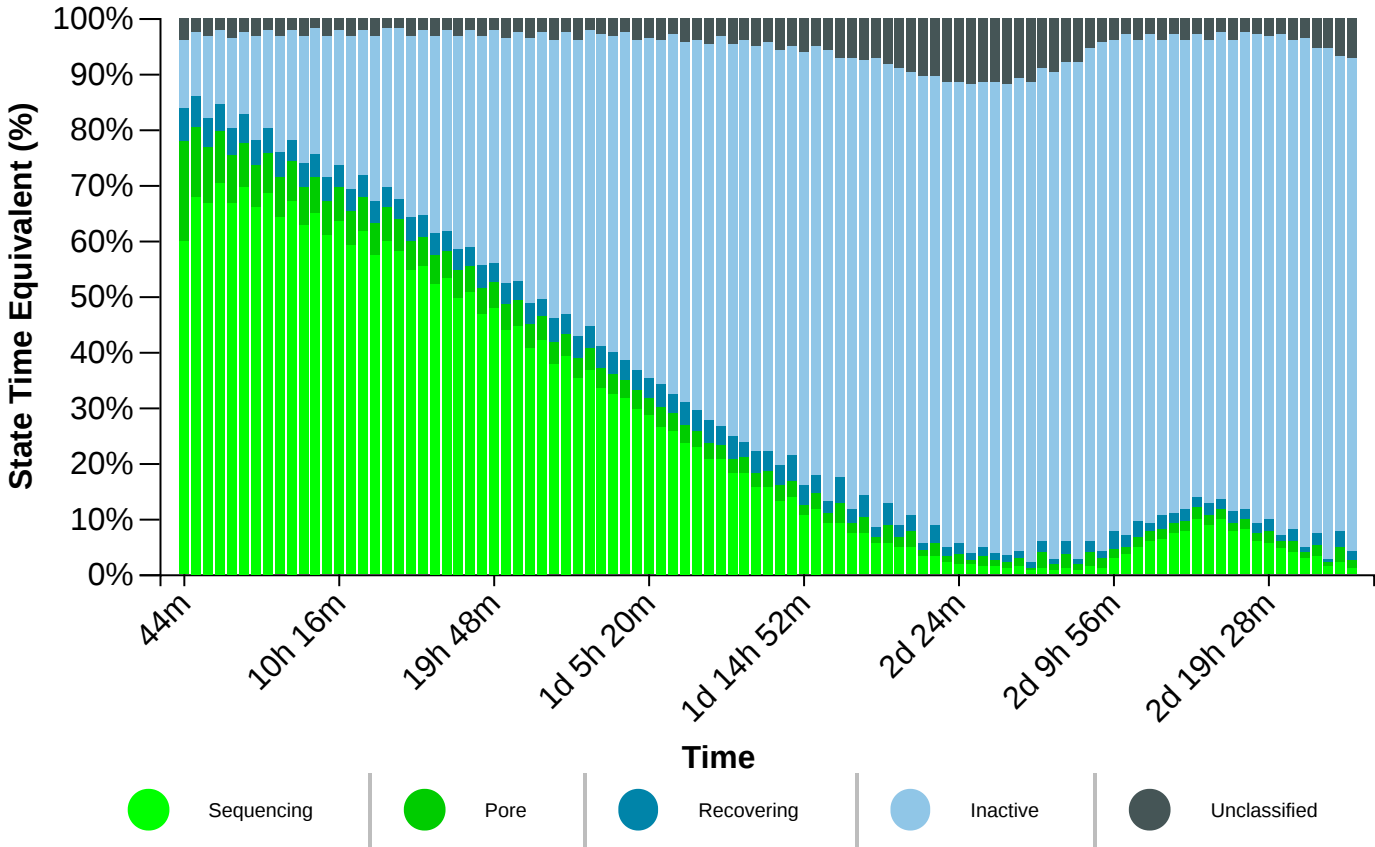


Duty time Categorised

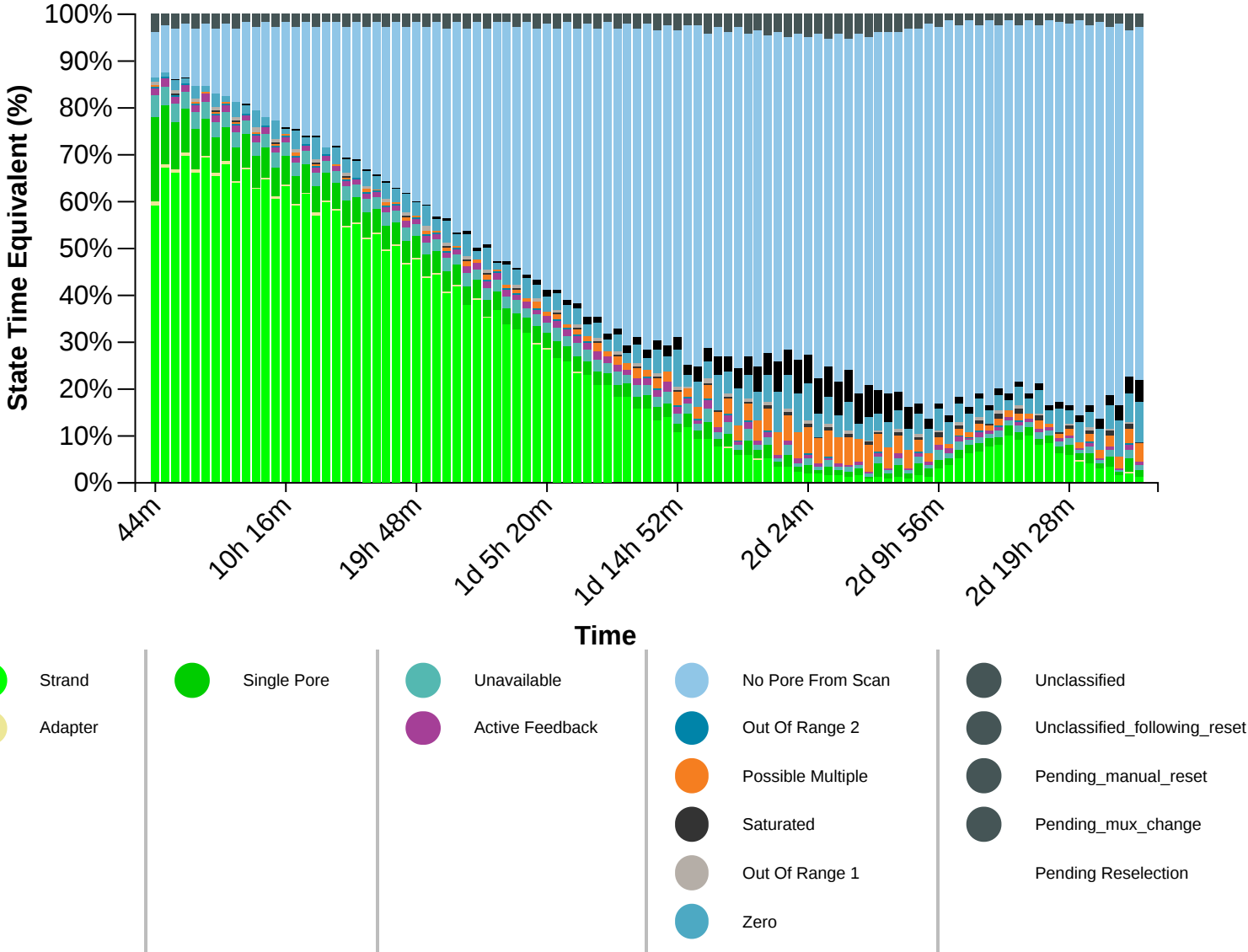


Mux Scan Grouped

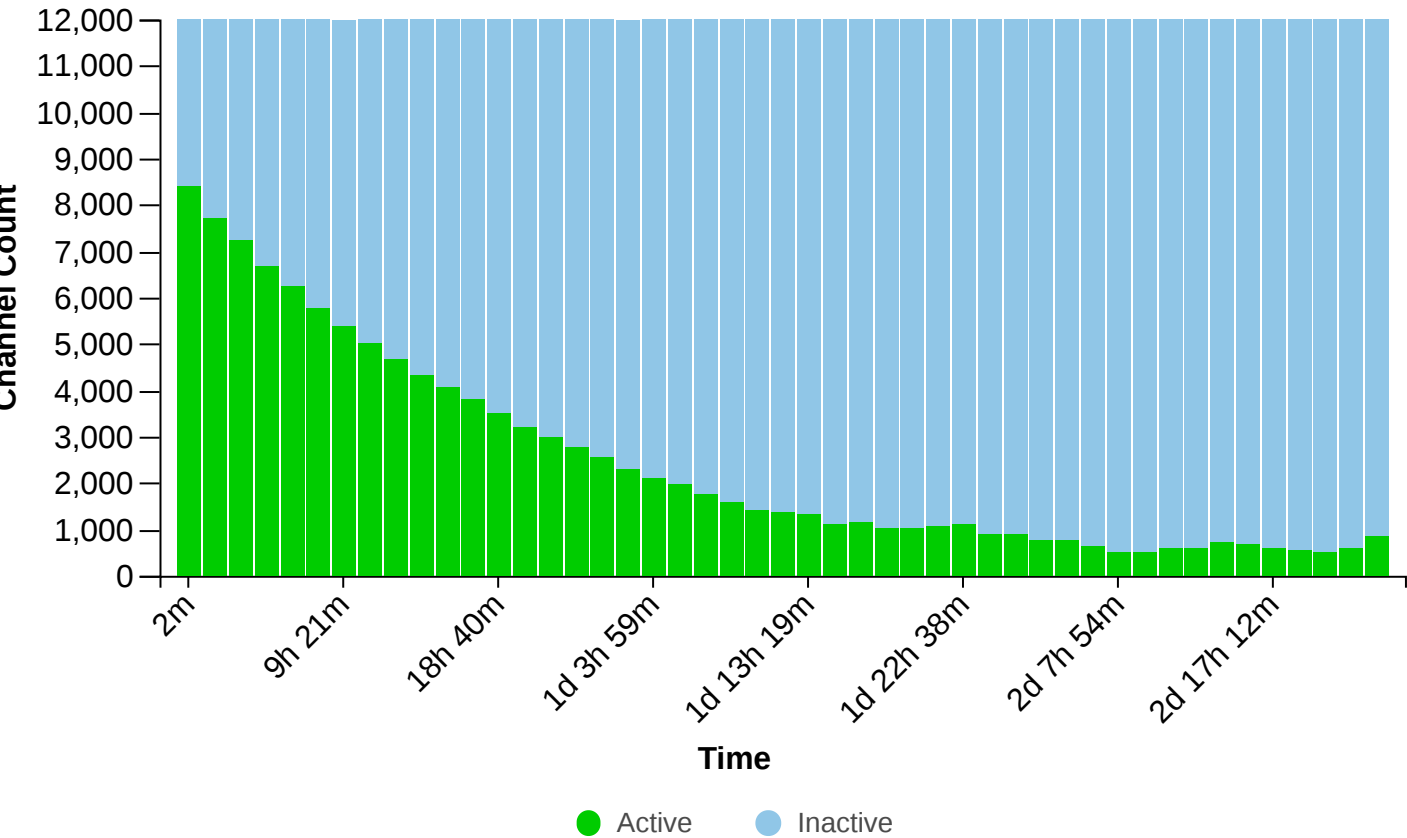


Mux Scan Categorised

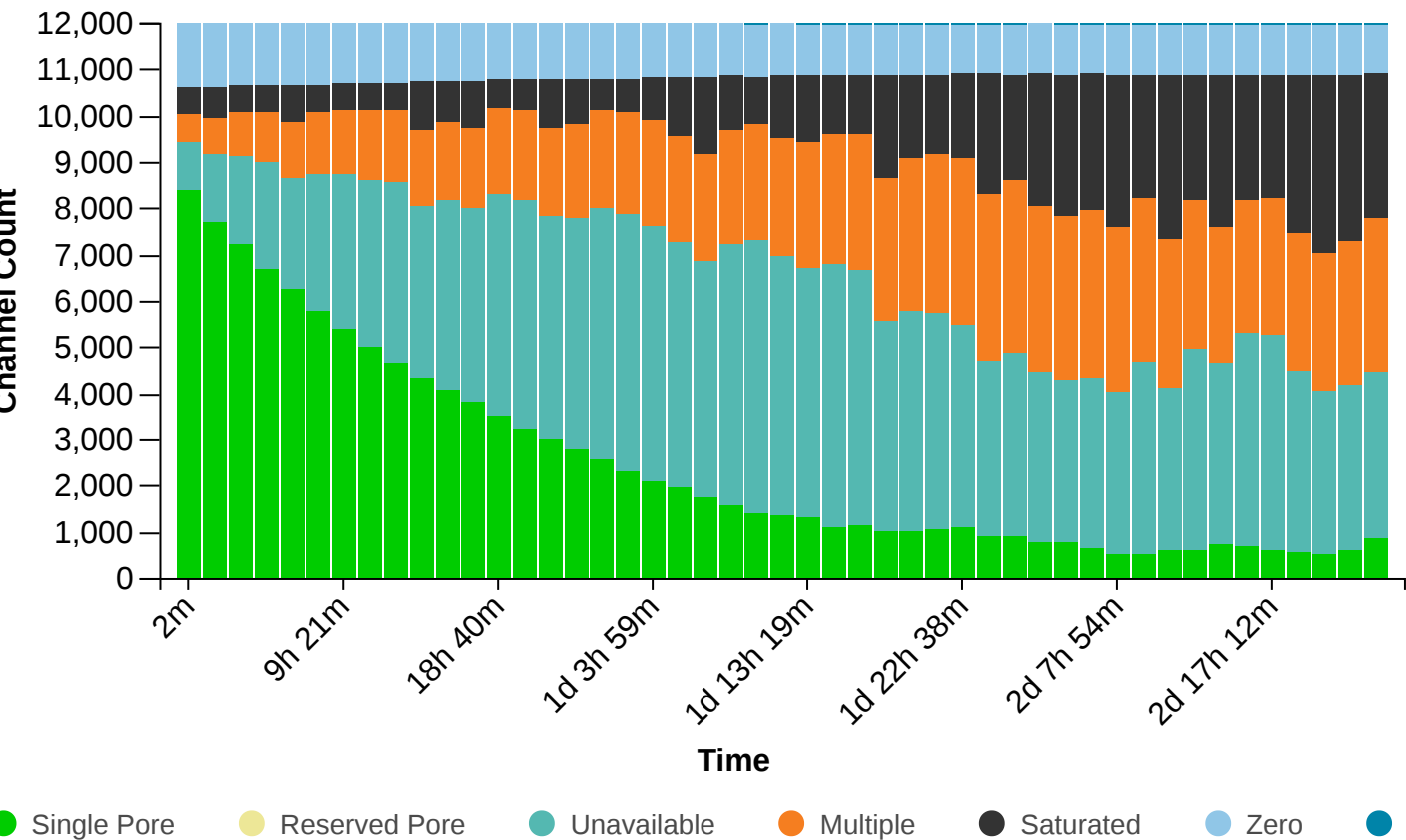


Temperature History

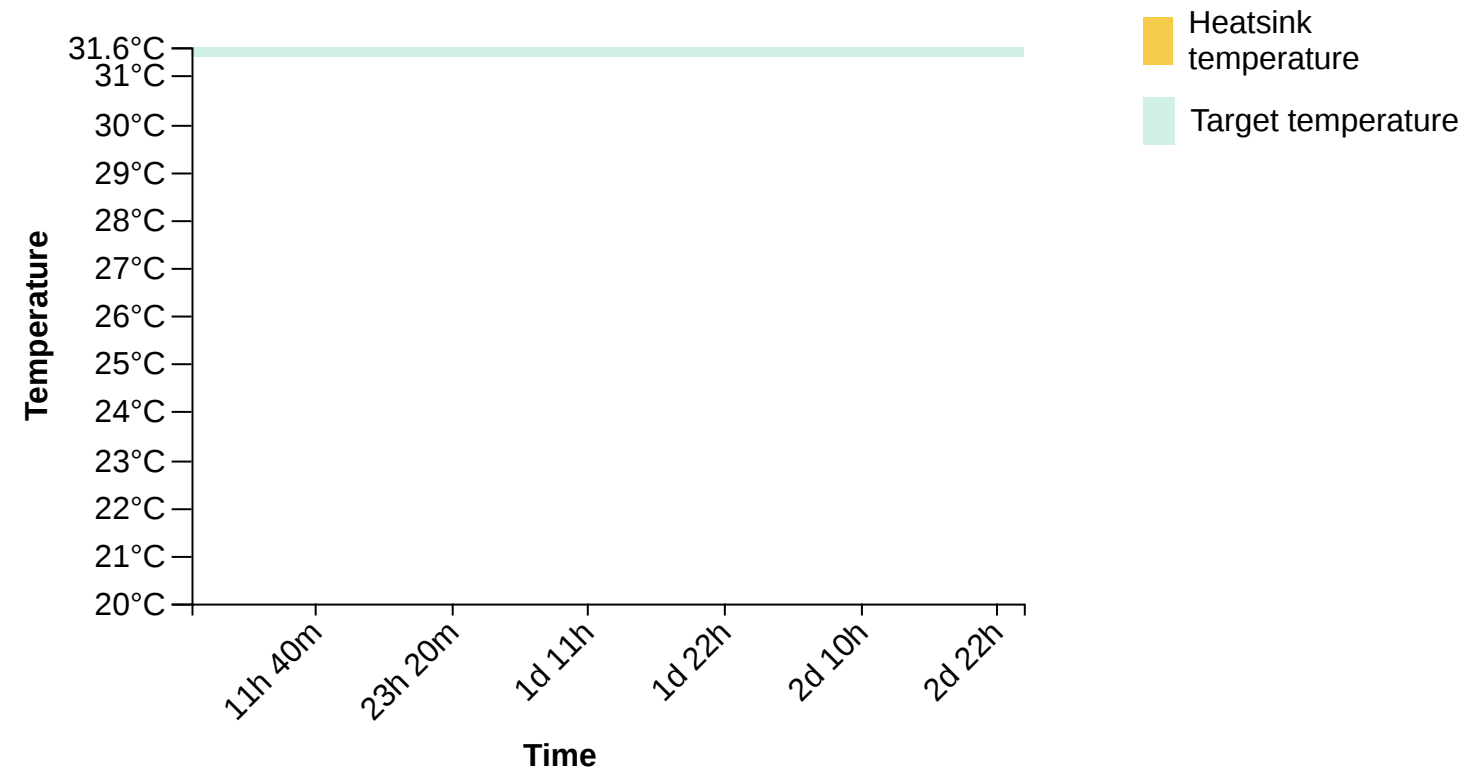


Bias Voltage History

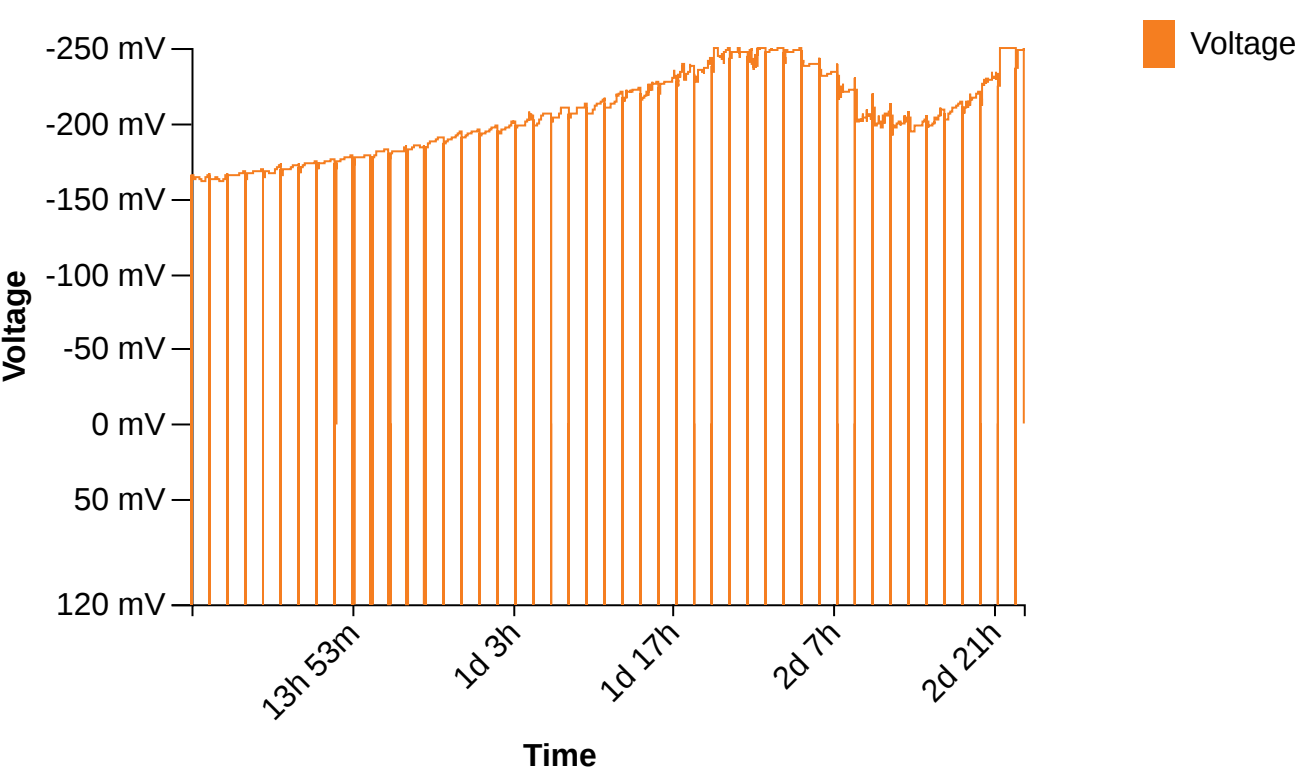


Translocation Speed

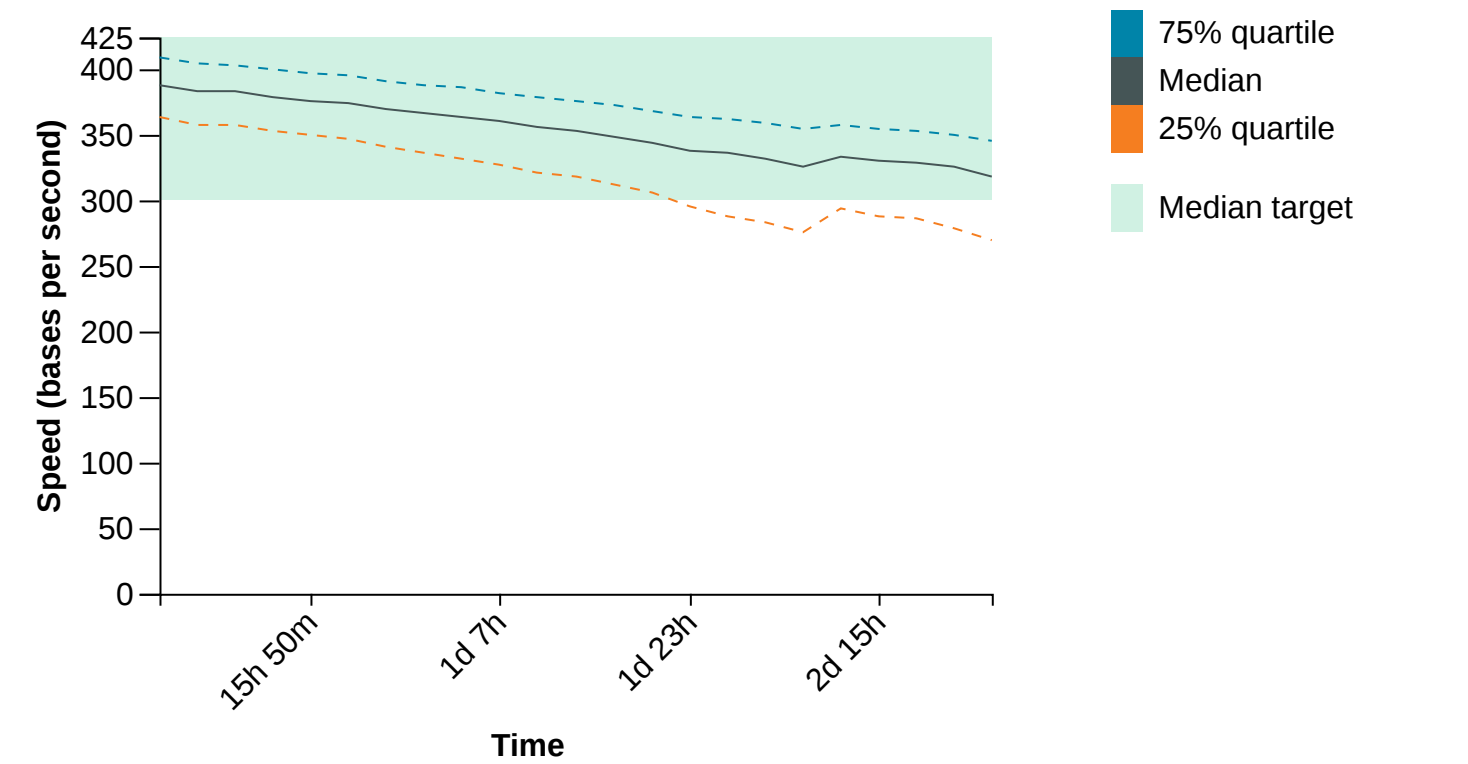

QScore

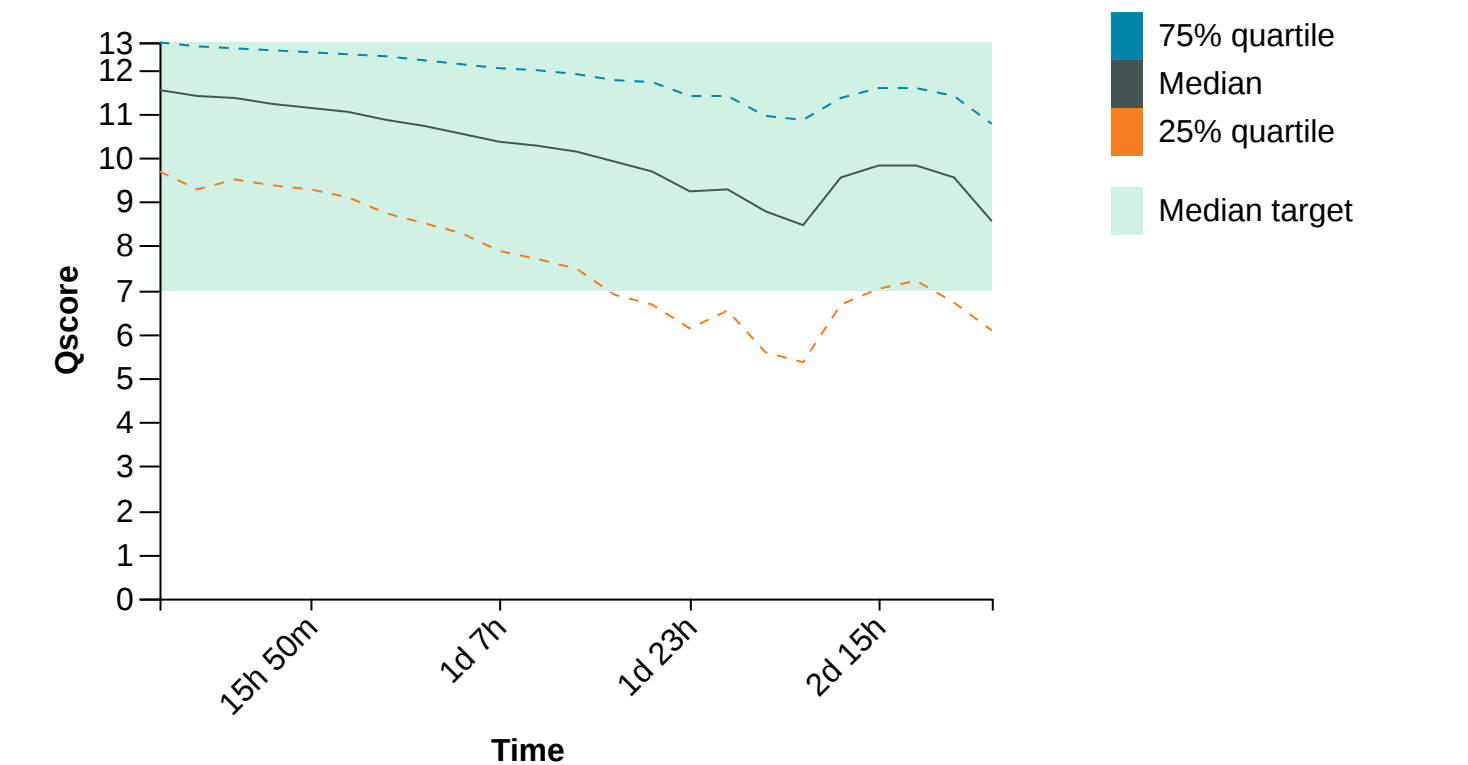

Disk Write Performance

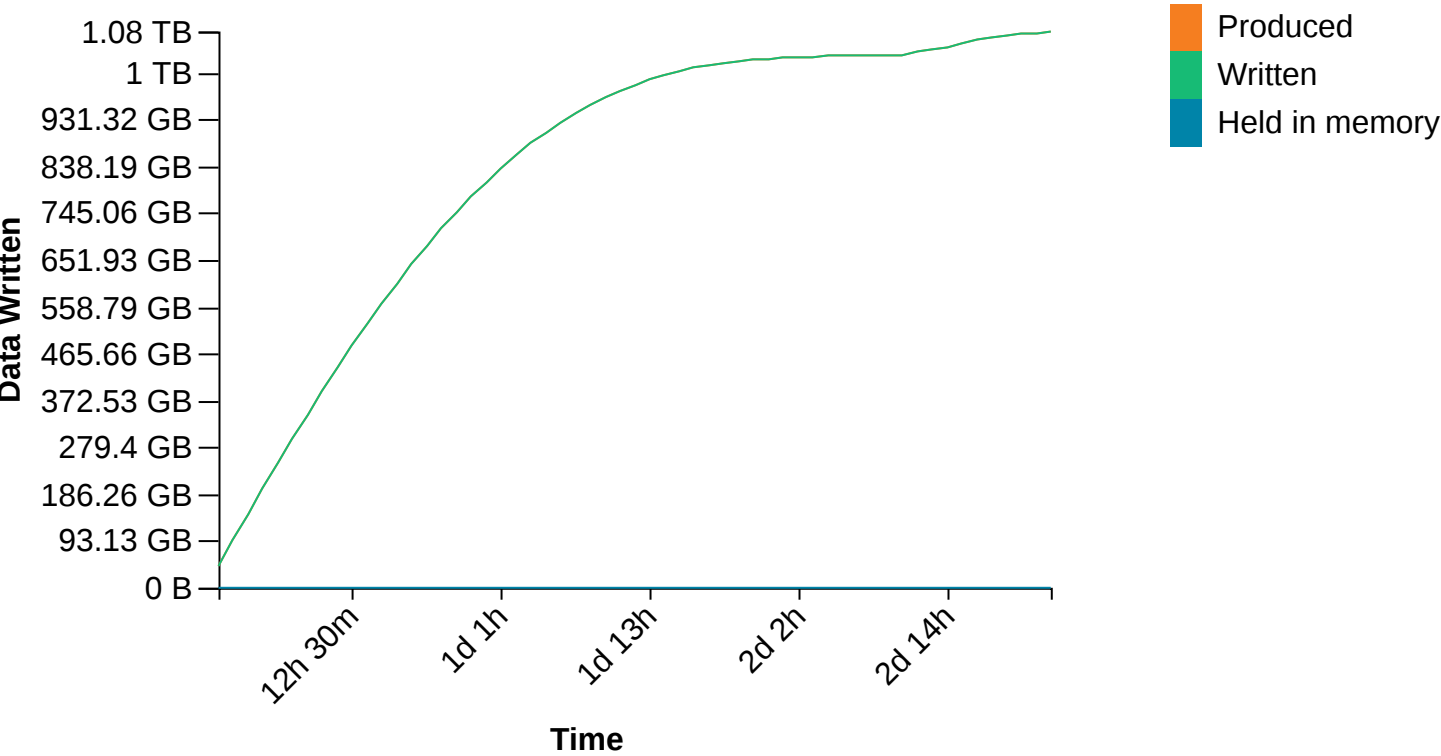

Run Debug Messages

- The sequencing run has finished, but basecalling may continue January 26, 16:06
- Flow cell PAE48907 has 836 pores available for sequencing. Starting sequencing with 732 pores January 26, 15:30
- Performing Mux Scan January 26, 15:27
- Flow cell PAE48907 has 602 pores available for sequencing. Starting sequencing with 556 pores January 26, 13:57
- Performing Mux Scan January 26, 13:54
- Flow cell PAE48907 has 505 pores available for sequencing. Starting sequencing with 459 pores January 26, 12:24
- Performing Mux Scan January 26, 12:21
- Flow cell PAE48907 has 538 pores available for sequencing. Starting sequencing with 470 pores January 26, 10:51
- Performing Mux Scan January 26, 10:48
- Flow cell PAE48907 has 608 pores available for sequencing. Starting sequencing with 535 pores January 26, 09:18
- Performing Mux Scan January 26, 09:15
- Flow cell PAE48907 has 692 pores available for sequencing. Starting sequencing with 611 pores January 26, 07:45
- Performing Mux Scan January 26, 07:42
- Flow cell PAE48907 has 709 pores available for sequencing. Starting sequencing with 614 pores January 26, 06:11
- Performing Mux Scan January 26, 06:08
- Flow cell PAE48907 has 621 pores available for sequencing. Starting sequencing with 545 pores January 26, 04:38
- Performing Mux Scan January 26, 04:35
- Flow cell PAE48907 has 586 pores available for sequencing. Starting sequencing with 529 pores January 26, 03:05
- Performing Mux Scan January 26, 03:02
- Flow cell PAE48907 has 511 pores available for sequencing. Starting sequencing with 473 pores January 26, 01:32
- Performing Mux Scan January 26, 01:29
- Flow cell PAE48907 has 503 pores available for sequencing. Starting sequencing with 467 pores January 25, 23:59
- Performing Mux Scan January 25, 23:57
- Flow cell PAE48907 has 654 pores available for sequencing. Starting sequencing with 580 pores January 25, 22:27
- Performing Mux Scan January 25, 22:24
- Flow cell PAE48907 has 766 pores available for sequencing. Starting sequencing with 677 pores January 25, 20:54
- Performing Mux Scan January 25, 20:52
- Flow cell PAE48907 has 791 pores available for sequencing. Starting sequencing with 691 pores January 25, 19:22
- Performing Mux Scan January 25, 19:19
- Flow cell PAE48907 has 910 pores available for sequencing. Starting sequencing with 773 pores January 25, 17:49
- Performing Mux Scan January 25, 17:47
- Flow cell PAE48907 has 915 pores available for sequencing. Starting sequencing with 795 pores January 25, 16:17

- Performing Mux Scan January 25, 16:14
- Flow cell PAE48907 has 1095 pores available for sequencing. Starting sequencing with 910 pores January 25, 14:44
- Performing Mux Scan January 25, 14:40
- Flow cell PAE48907 has 1076 pores available for sequencing. Starting sequencing with 887 pores January 25, 13:10
- Performing Mux Scan January 25, 13:07
- Flow cell PAE48907 has 1045 pores available for sequencing. Starting sequencing with 846 pores January 25, 11:37
- Performing Mux Scan January 25, 11:34
- Flow cell PAE48907 has 1041 pores available for sequencing. Starting sequencing with 815 pores January 25, 10:04
- Performing Mux Scan January 25, 10:01
- Flow cell PAE48907 has 1153 pores available for sequencing. Starting sequencing with 892 pores January 25, 08:31
- Performing Mux Scan January 25, 08:28
- Flow cell PAE48907 has 1121 pores available for sequencing. Starting sequencing with 818 pores January 25, 06:58
- Performing Mux Scan January 25, 06:55
- Flow cell PAE48907 has 1341 pores available for sequencing. Starting sequencing with 954 pores January 25, 05:24
- Performing Mux Scan January 25, 05:21
- Flow cell PAE48907 has 1362 pores available for sequencing. Starting sequencing with 918 pores January 25, 03:51
- Performing Mux Scan January 25, 03:48
- Flow cell PAE48907 has 1423 pores available for sequencing. Starting sequencing with 941 pores January 25, 02:18
- Performing Mux Scan January 25, 02:15
- Flow cell PAE48907 has 1575 pores available for sequencing. Starting sequencing with 1013 pores January 25, 00:45
- Performing Mux Scan January 25, 00:42
- Flow cell PAE48907 has 1757 pores available for sequencing. Starting sequencing with 1125 pores January 24, 23:12
- Performing Mux Scan January 24, 23:08
- Flow cell PAE48907 has 1990 pores available for sequencing. Starting sequencing with 1224 pores January 24, 21:38
- Performing Mux Scan January 24, 21:35
- Flow cell PAE48907 has 2119 pores available for sequencing. Starting sequencing with 1297 pores January 24, 20:05
- Performing Mux Scan January 24, 20:02
- Flow cell PAE48907 has 2326 pores available for sequencing. Starting sequencing with 1382 pores January 24, 18:32
- Performing Mux Scan January 24, 18:29
- Flow cell PAE48907 has 2557 pores available for sequencing. Starting sequencing with 1475 pores January 24, 16:59
- Performing Mux Scan January 24, 16:56
- Flow cell PAE48907 has 2765 pores available for sequencing. Starting sequencing with 1553 pores January 24, 15:26
- Performing Mux Scan January 24, 15:22
- Flow cell PAE48907 has 3011 pores available for sequencing. Starting sequencing with 1651 pores January 24, 13:52
- Performing Mux Scan January 24, 13:49
- Flow cell PAE48907 has 3230 pores available for sequencing. Starting sequencing with 1758

- pores January 24, 12:19
- Performing Mux Scan January 24, 12:16
- Flow cell PAE48907 has 3530 pores available for sequencing. Starting sequencing with 1849 pores January 24, 10:46
- Performing Mux Scan January 24, 10:43
- Flow cell PAE48907 has 3793 pores available for sequencing. Starting sequencing with 1940 pores January 24, 09:13
- Performing Mux Scan January 24, 09:10
- Flow cell PAE48907 has 4064 pores available for sequencing. Starting sequencing with 2028 pores January 24, 07:39
- Performing Mux Scan January 24, 07:36
- Flow cell PAE48907 has 4316 pores available for sequencing. Starting sequencing with 2119 pores January 24, 06:06
- Performing Mux Scan January 24, 06:03
- Flow cell PAE48907 has 4692 pores available for sequencing. Starting sequencing with 2200 pores January 24, 04:33
- Performing Mux Scan January 24, 04:30
- Flow cell PAE48907 has 5016 pores available for sequencing. Starting sequencing with 2267 pores January 24, 03:00
- Performing Mux Scan January 24, 02:57
- Flow cell PAE48907 has 5396 pores available for sequencing. Starting sequencing with 2329 pores January 24, 01:27
- Performing Mux Scan January 24, 01:24
- Flow cell PAE48907 has 5783 pores available for sequencing. Starting sequencing with 2395 pores January 23, 23:53
- Performing Mux Scan January 23, 23:50
- Flow cell PAE48907 has 6264 pores available for sequencing. Starting sequencing with 2478 pores January 23, 22:20
- Performing Mux Scan January 23, 22:17
- Flow cell PAE48907 has 6703 pores available for sequencing. Starting sequencing with 2539 pores January 23, 20:47
- Performing Mux Scan January 23, 20:44
- Flow cell PAE48907 has 7240 pores available for sequencing. Starting sequencing with 2608 pores January 23, 19:14
- Performing Mux Scan January 23, 19:11
- Flow cell PAE48907 has 7734 pores available for sequencing. Starting sequencing with 2655 pores January 23, 17:41
- Performing Mux Scan January 23, 17:38
- Flow cell PAE48907 has 8392 pores available for sequencing. Starting sequencing with 2702 pores January 23, 16:07
- Performing Mux Scan January 23, 16:05
- Starting sequencing procedure January 23, 16:05
- 180 seconds have elapsed. Experiment commencing. January 23, 16:04
- Waiting up to 180 seconds for temperature to stabilise at 31.5°C January 23, 16:01
- Disk /data has 16165 GB space remaining January 23, 16:01